



Vet-Tables

By Dispomed

Vet-Tables Scissor Table



Multi task table

Whenever you need a stretcher or extra help to lift a heavy dog or perform surgery, the scissor table is the perfect table to help you accomplish your daily tasks.

Essential For Any Veterinary Clinic

The exceptional lifting capacity of the Vet-Tables scissor table protects your patients and staff from injuries caused by the handling of heavy animals.

- Lift capacity 300 pounds (136 Kg)
- Canadian made
- Rapid rail tie down system lets you safely position animals
- The mobile parts of the table are on bearings
- Every type of table top is made of 16 gauge, 304 grade stainless steel.
- Minimum height 11" (27.9 cm)
- Maximum height 40" (101.6 cm)
- Slide cleaners prevent the accumulation of hair in the lifting rails
- Full aluminium structure for a greater robustness
- The silent motor stops automatically when the maximum height is reached
- 5 year warranty

Vet-Tables Scissor Table



You need your table to be **customized to your needs**, and for this reason Dispomed can propose customized options on anything from the base to the top and anything in between.

1. Table Bases Options

Mobile Base



The mobile base includes 4" wheels with brakes. The mobile base also comes with a battery in order to use it when disconnected. the battery, when fully charged, has power for 30 lifts.

Fixed Base

The fixed base is perfect for an examination room, giving a secure way to position heavy dogs on a table with a very stable base.



2. Table Top Options



Exam top

Standard exam top is perfect for an examination room. This top is available in 4 and 5 feet.

Surgery top

Indented contour with a hole and a hook at one extremity permitting the accumulation of liquids during surgery. The surgery top includes 4 paw attachments that can be moved on the 2 rails on each side of the table. Available in 4 and 5 feet.



Vet-Tables Scissor Table



V top

The V top gives you an adjustment with 2 independent handles that allow you to adjust the panels independently. There is a drawer between the 2 panels to accumulate liquids and waste from surgery. The V top includes 4 paw attachments that can be moved on the 2 rails on each side of the table. Available in 4 and 5 feet.



Dental top

The dental top is the best top to use when performing a dentistry procedure. The easy to remove grill allows the liquid to evacuate into the basin. The drain valve is under the basin and is easy to access, which allows for fast and easy cleaning of the dental table. The dental top includes 4 moveable paw attachments placed on the rails on each side of the table. The foot pedal is attached to the front of the table for easy access when performing dentistry procedures.



3. More Options

Once you have decided which base and top you want, there are still a couple of options you can add to your customized table.

Tilt

All of our tops are available with the electric tilt option permitting you to have the perfect angle to work. The tilt option comes with a double electrical pedal so your hands remain free.



Heated

If you want extra help to avoid heat loss, add a heated top. This is the easiest way to keep your patient warm. The top is heated with a water circuit, so no risk of electric hazards as there is no temperature difference from one extremity to the other. This option comes with the HTP 1500 water pump. This option is not available on the V top.




Vet-Tables Scissors Table



Color option

You can match your customized table with your clinic, by choosing the perfect color and creating the perfect combo with our rock gray standard color.

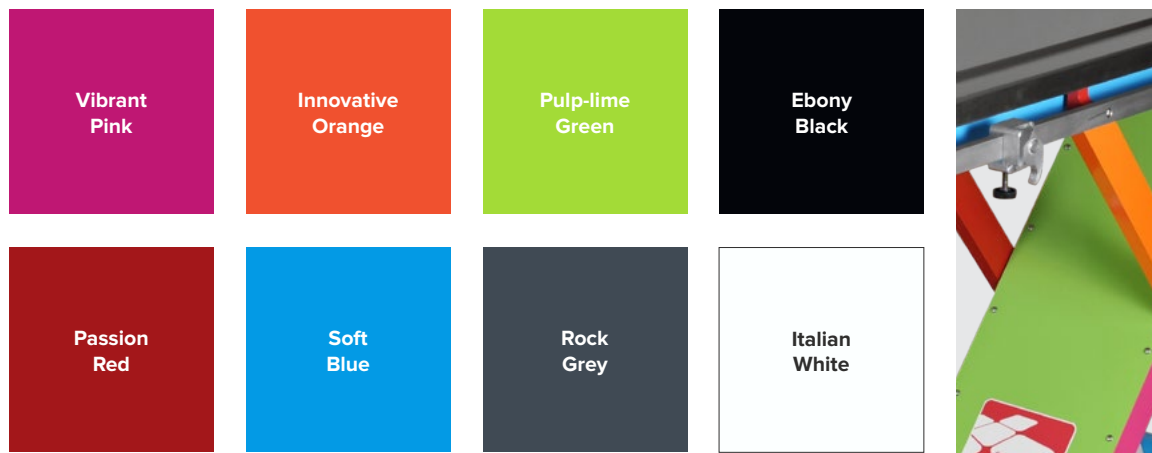

Available Configurations

Fixed

4' (1.21m) exam top : # 285-0200-4EXX
4' (1.21m) surgery top heated, tilt : # 285-0200-4SHT
4' (1.21m) surgery top heated : # 285-0200-4SHX
4' (1.21m) surgery top tilt : # 285-0200-4SXT
4' (1.21m) surgery top : # 285-0200-4SXX
4' (1.21m) "V" top tilt : # 285-0200-4VXT
4' (1.21m) "V" top : # 285-0200-4VXX
5' (1.52m) exam top : # 285-0200-5EXX
5' (1.52m) surgery top heated tilt : # 285-0200-5SHT
5' (1.52m) surgery top heated : # 285-0200-5SHX
5' (1.52m) surgery top tilt : # 285-0200-5SXT
5' (1.52m) surgery top : # 285-0200-5SXX
5' (1.52m) "V" top tilt : # 285-0200-5VXT
5' (1.52m) "V" top : # 285-0200-5VXX
5' (1.52m) dental top : # 285-0200-5DXX
5' (1.52m) dental top, tilt : # 285-0200-5DXT
5' (1.52m) dental top, tilt, heated : # 285-0200-5DHT
5' (1.52m) dental top, heated : # 285-0200-5DHX

Mobile

4' (1.21m) exam top : # 285-0210-4EXX
4' (1.21m) surgery top, heated, tilt : # 285-0210-4SHT
4' (1.21m) surgery top, heated : # 285-0210-4SHX
4' (1.21m) surgery top, tilt : # 285-0210-4SXT
4' (1.21m) surgery top : # 285-0210-4SXX
4' (1.21m) "V" top, tilt : # 285-0210-4VXT
4' (1.21m) "V" top : # 285-0210-4VXX
5' (1.52m) exam top : # 285-0210-5EXX
5' (1.52m) surgery top, heated, tilt : # 285-0210-5SHT
5' (1.52m) surgery top, heated : # 285-0210-5SHX
5' (1.52m) surgery top, tilt : # 285-0210-5SXT
5' (1.52m) surgery top : # 285-0210-5SXX
5' (1.52m) "V" top, tilt : # 285-0210-5VXT
5' (1.52m) "V" top : # 285-0210-5VXX
5' (1.52m) dental top : # 285-0210-5DXX
5' (1.52m) dental top, tilt : # 285-0210-5DXT
5' (1.52m) dental top, tilt, heated : # 285-0210-5DHT
5' (1.52m) dental top, heated : # 285-0210-5DHX

Vet-Tables Electric Surgery Table



Small footprint

If you are looking for a table with the most room for your legs, you certainly have to look at our electric column table.

The Vet-Tables electric column table has a large possibility of configurations. This table is manufactured in Canada and designed exclusively for veterinary use. The electric column table offers stability with a minimum footprint.

- 5 year warranty
- Manufactured in Canada
- Lift capacity of 300 pounds (136 Kg)
- Every type of table top is made of 16 gauge, 304 grade stainless steel.
- Rapid rail tie down system lets you safely position animals
- Maximum height of 43,5 in.
- Minimum height of 32,25 in.
- Silent motor
- Electric components are water resistant
- Manual tilt
- 360 degree turntable
- Table top locking system for greater stability
- Weight: 305 lbs (138.3 Kg).

Vet-Tables Electric Surgery Table



1. Table Top Options

Surgery top

Indented contour with a hole and a hook at one extremity permitting the accumulation of liquids during surgery. The surgery top includes 4 paw attachments that can be moved on the 2 rails on each side of the table. Available in 4 and 5 feet.



V top

The V top gives you an adjustment with 2 independent handles that allow you to adjust the panels independently. There is a drawer between the 2 panels to accumulate liquids and waste from surgery. The V top includes 4 paw attachments that can be moved on the 2 rails on each side of the table. Available in 4 and 5 feet.



Vet-Tables Electric Surgery Table



Dental Top Options

The dental top is the best top to use when performing a dentistry procedure. The easy to remove grill allows the liquid to evacuate into the basin. The drain valve is under the basin and is easy to access, which allows for fast and easy cleaning of the dental table. The dental top includes 4 moveable paw attachments placed on the rails on each side of the table.



More Options

Once you have decided which base and top you want, there are still a couple of options you can add to your customized table.

Heated top

If you want extra help to avoid the heat loss, add a heated top. This is the easiest way to keep your patient warm. The top is heated with a water circuit, so no risk of electric hazards as there is no temperature difference from one extremity to the other. This option comes with the HTP 1500 water pump. This option is not available on the V top.

Vet-Tables Electric Surgery Table



Available Configurations

4' (1.21m) surgery top, heated : # 285-0400-4SHT

4' (1.21m) surgery top : # 285-0400-4SXT

4' (1.21m) V top : # 285-0400-4VXT

5' (1.52m) surgery top, heated : # 285-0400-5SHT

5' (1.52m) surgery top : # 285-0400-5SXT

5' (1.52m) V top : # 285-0210-4SHX

5' (1.52m) Dental top : # 285-0400-5DXT

5' (1.52m) Dental top, heated : # 285-0400-5DHT

I.V. pole for Vet-Tables Electric and scissor tables

The height for our optional I.V. pole for our scissor and electric tables is adjustable. It is practical to hang your infusion pump or any other accessories you want to have close to your patient.





Extra workspace

If you need extra space for your very busy day at a low price, this table is for you.

H base exam table

- One piece welded.
- Type 304 stainless steel construction.
- Flat top
- Very stable

1. Length

Choose between a 4 or 5 foot long top.

Vet-Tables Exam Table



2. Base



Fixed

The fixed exam table is perfect for a designated space



Mobile

The mobile base includes 4" wheels with brakes.

Available Configurations

Fixed

4' (1.21m): # 285-0548-000

5' (1.52m): # 285-0560-000

Mobile

4' (1.21m): # 285-0545-000

5' (1.52m): # 285-0565-000

Vet-Tables Mobile Exam and Surgery Table



Lack of space

If you lack space, our mobile anesthesia table trolley is the perfect combo for you. You can use it as a stretcher as well as for surgery purposes.

Mobil anesthesia table trolley

- Available with or without anesthesia machine.
- 4" locking castors.
- Built-in bumpers.
- 2 shelves.
- Liquid drain & pail hook.
- Drawer assembly shown (optional).

Contact your representative for specific configurations.

The Vet-Tables are proudly Made in Canada





Made by **DISPOMED**



745 Nazaire-Laurin, Joliette (Qc) Canada J6E 0L6
Tel: 1 800 363-1746 - Fax: 450 759-8181
info@dispomed.com - www.dispomed.com