

# LifeWindow™ One

Veterinary Multi-Parameter Physiologic Monitor

## Monitor Confidently



### Gold standard veterinary measurements

- 3/5 Lead Digicare™ VetECG™
- Masimo® SET® Rainbow® oximetry with PVI®
- Capnography with optional Anesthetic Agents
- SunTech® VetBP® Blood Pressure
- Temperature
- Up to 8 hours battery life

*Driven by Innovation ...*  
*... in Veterinary Medicine*  
[www.digi-vet.com](http://www.digi-vet.com)



# LifeWindow™ One **Veterinary Multi-Parameter Physiologic Monitor**

Digicare's new LifeWindow™ One is a portable veterinary specific physiological monitor with distinctive features and benefits. Engineered with premium veterinary specific modules including an easily operated touch screen user interface. This monitor is in a class of its own, compact and weighing only 3 Lbs. Made with convenience in mind, the LifeWindow™ One can be easily mounted on a rolling stand or IV pole. The perfect combination of size and functionality to complement any hospital (clinic) setting.



Our most popular models are below. For a complete list of models please visit our website.

<b>LifeWindow One</b>					
	Code	Picture	Measurements	Modules	Description
Single (Includes Acc. for first param)	LW11		NIBP PnP #1 PnP #2	Suntech PnP Connector for SpO2 PnP Connector for CO2	NIBP only with connections for CO2 and SpO2
	LW21		SpO2 NIBP PnP	Masimo SpO2 Internal Suntech PnP Connector for CO2	Masimo SpO2, Suntech NIBP and connection for CO2
Dual (Includes accessories for first two parameter)	LW23		CO2 SpO2	Respironics (Sidestream Internal) Masimo SpO2 Internal	Respironics sidestream CO2 and Masimo SpO2
	LW26		SpO2 Multi-Gas	Masimo SpO2 Internal Masimo Gas (Sidestream Internal AX+)	Masimo SpO2 and Sidestream gas detection
	LW27		SpO2 Multi-Gas	Masimo SpO2 Internal Masimo Gas (Mainstream AX+)	Masimo SpO2 and Mainstream gas detection
Triple (Includes All Accessories)	LW31		ECG SpO2 NIBP	Digicare Masimo SpO2 Internal Suntech	ECG, Masimo SpO2, and Suntech NIBP
	LW33		NIBP CO2 SpO2	Suntech Respironics (Sidestream Internal) Masimo SpO2 External Module	Suntech NIBP, Respironics sidestream internal CO2, Masimo SpO2 (Plug and play)

**\* Notes \***


1. Temperature is included in all models, accessory sold separately.
2. PnP modules - Monitor is capable of measuring the parameter, accessory is sold separately.
3. PnP SpO2 or CO2 - Monitor is pre-configured for both CO2 and SpO2, decision does not need to be made at time of sale.
4. Masimo SpO2 Internal - SpO2 Module is installed inside of the monitor, this is a Rainbow module and can accept PVI + other Rainbow Parameters.
5. Masimo SpO2 External - SpO2 Module is installed outside of the monitor, this module cannot have PVI or any other Rainbow Parameter.
6. Other measurement module options available upon request, check with Digicare for availability.

**General Specifications**

**Display Type** - 480x272 4.8" LED backlit TFT screen  
**Mains Voltage** - Automatic Selection 100Vac to 240Vac (50/60 Hz)  
**Battery** - Internal Li-ion rechargeable battery

**Relative Humidity** - 30 to 70% (non condensing)  
**Dimensions** - W 9.14"(232) H 3.46"(88) D 5.67"(144) (mm)  
**Weight** - 3 Lbs (1.4 Kg)  
**Safety** - Meets IEC601-1-1 and IEC 601-1-2

**Driven by Innovation in Veterinary Medicine**

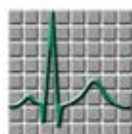

**DIGICARE**  
Animal Health

Digicare Animal Health  
 107 Commerce Road  
 Boynton Beach, Florida, 33426  
 USA  
 Phone: 561-689-0408  
 Fax: 561-689-0021

Distributed by:



For more information go to:  
**WWW.DIGI-VET.COM**  
 email:sales@digicarebiomedical.com

**844-4DIGICARE**